

TEL (02) 533-6720
서울시 동작구 신대방1가길 38 (신대방 719 동작상떼빌) 106동 209호

DTMc300 Double Monochromator

www.wonwoosystem.co.kr

DTMc300 Monochromator

For applications where the scattered light performance of a single monochromator is not sufficient, the excellent optical performance as the TMc300 may still be experienced with two of such used in tandem constituting thereby the DTMc300 double monochromator.


The DTMc300 maintains the flexibility of configuration and ability to cover a wide spectral range of the TMc300, within a single scan.

Compatible with the full range of sources, detectors and accessories of the Bentham portfolio, in the additive dispersion mode, the DTMc300 monochromator is at the heart of many of our light measurement solutions, particularly those incorporating a UV component, whilst the subtractive dispersion mode is sought by high level laboratories seeking minimum measurement uncertainty due to the resulting lack of dispersion at the exit slit.

Key features:


- Double monochromator composed of two 300mm focal length monochromators in Czerny-Turner configuration arranged for either additive or subtractive dispersion
- Each single monochromator accommodates up to three diffraction gratings permitting single-scan measurement over wide spectral range included in range 200nm-50µm
- Up to two entrance ports, one exit port on first monochromator (single mode) and up to two exits on second monochromator (double mode) with automated mirror-based port selection for use of multiple input/ output/ detector configurations
- Fully automated through USB 2.0 interface

Configuration (typical)

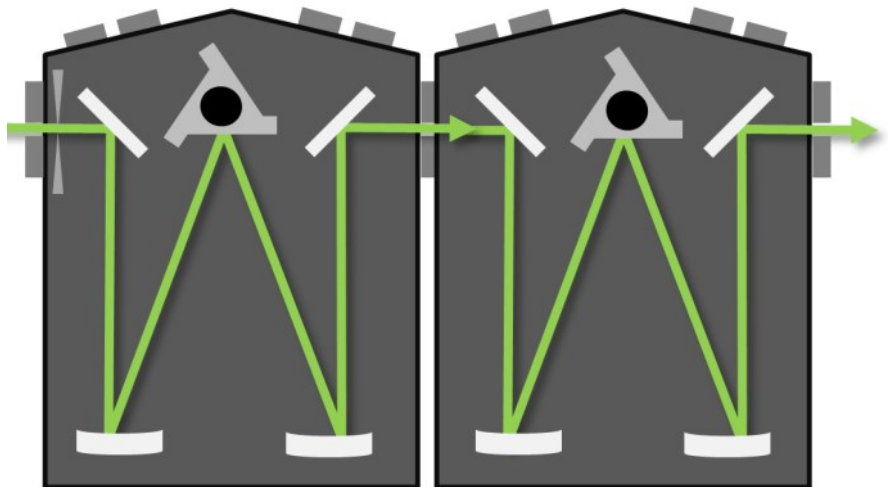

Standard configuration A

Alternative configurations B & C:

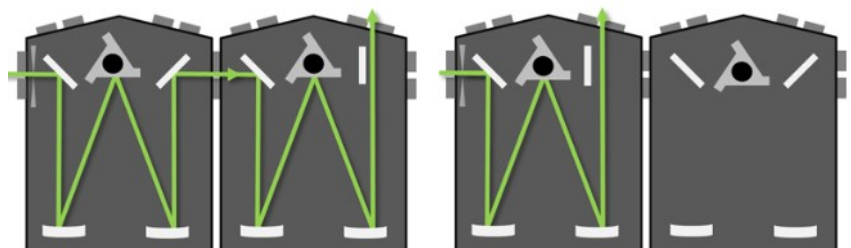


Specification


Optical Layout	
Configuration	Pair symmetrical Czerny Turner monochromators arranged for either additive or subtractive dispersion
Focal length	Total focal length: 600mm (additive) Individual Units: 300mm
Grating mount	Triple grating turret
Number of gratings	1-3
Grating size	68 x 84mm
Aperture ratio	f/4.1 (at all grating angles)
Number of entrance/exit ports	1 or 2 in both monochromators
Opto-Mechanical	
Mechanical resolution drive	0.000072°/ step
Slit type	Fixed, micrometer or motorised variable
Slit width & height	10µm-10mm (W) x 20mm (H)
Filter size	23/25mm diameter
Optical Performance (quoted for 1200 g/mm gratings)	
Spectral range	Within range 200nm-50µm
Linear dispersion	1.35nm/mm (additive)
Wavelength accuracy	± 0.1nm
Wavelength reproducibility	± 0.025nm
Resolution full/reduced slit height	0.15/0.05nm
Stray light rejection at 2.5 FWHM	10 ⁻⁸
Wavelength acquisition speed	2000nm/s
Automation	
Computer interface	USB
Software control	BenWin+/SDK
Electrical/Mechanical	
Overall dimensions (mm)	634W x 360D x 211H
Optical height (mm)	115mm
Weight	28 kg
Power supply	Mains input 110/220V 50/60Hz
Options	
Diffraction gratings	See separate Diffraction Gratings datasheet
Order sorting filters	See separate Order Sorting Filters datasheet
Filter wheel	252-6 or 252-8 (6/8-position filter wheel,



DTMc300 Dimensions

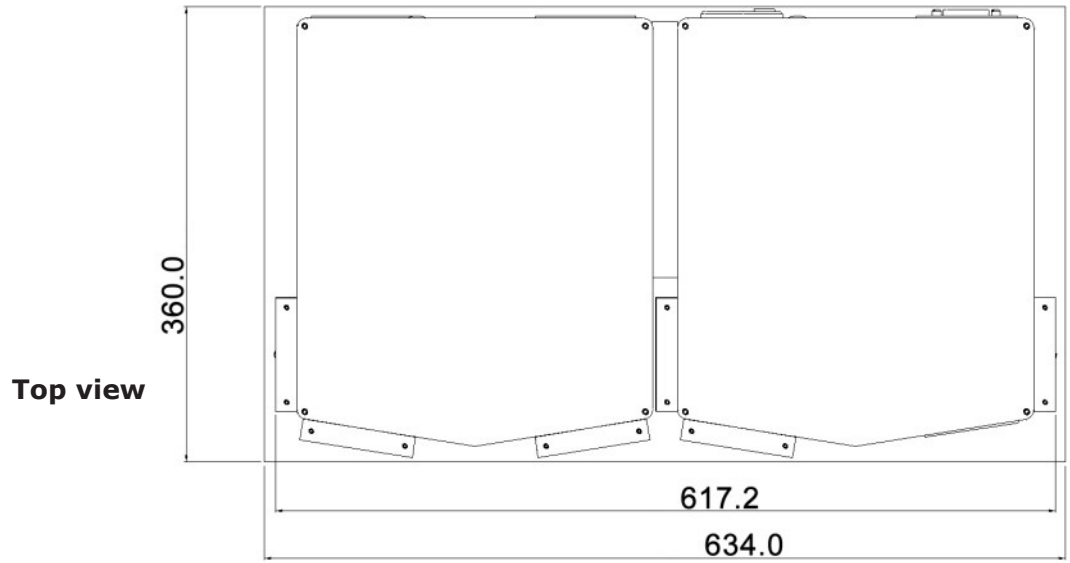

Side view

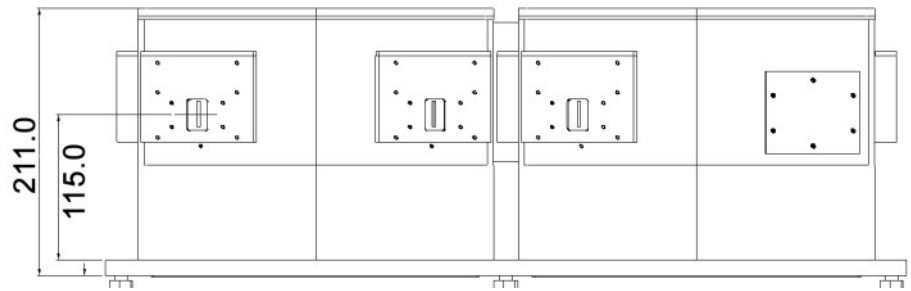

DTMc300V3.1